

Warranty Guidelines

- Visit <https://smartskn.com/product-registration> to access your warranty and register your device.
- The Muilli device has a **1-year warranty from the date of purchase** when used as directed
- The warranty does not cover damages from misuse, unauthorized repairs, or accidents, only manufacturing defects.
- Use only authentic parts and follow the product manual carefully.
- You may return defective devices with proof of purchase for a refund or replacement under Warranty Conditions.

Service Guidelines

Service is covered for issues that arise when using as directed for one full year after purchase, starting on the day of purchase.

For secondhand devices, the 1-year warranty begins at the time of the original purchase date.

Paid Service Guidelines

You will be charged for damages due to negligence or failure to follow instructions.

Damage from unauthorized repairs, improper electrical use, or natural disasters is not covered.

Important Safety Information for your Muilli device

Do not open, disassemble, or tamper with the device or its chargers to avoid malfunction and safety hazards.

Do not use water or liquid soaps or detergents to clean the device.

Do not use in saunas or extremely humid conditions.

Avoid leaving in direct sunlight for extended periods of time.

Use a charger with a 5.0V output voltage and at least 1A output current.

Keep the device away from fire to avoid risk of explosion.

Store the device out of the reach of children.

Health And Safety Precautions

The Muilli is not intended for medical use.

The Muilli cannot diagnose, treat, or prevent illnesses or diseases.

Please consult a healthcare professional for advice on skin conditions and treatment.

FCC Information to User

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Caution

Modification not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC Compliance Information

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this must accept any interference received, including interference that may cause undesired operation.

Product name	Muilli 2
Model	LC-200
Size	52.8 x 34.7 x 174
Weight	140g
Material	PC/multifunction PC
Rated voltage/power consumption	3.6V 300mA
Major specification	WIFI-IEEE 802.11 b/g/n (Duplex transmission)
Operating temperature	-20°C ~ 50°C
Charging specification	5V/1A
Battery	3.6V/3250mA Lithium-ion
Time needed for full charging	Charging time 5 hours Usage time 8 hours
Manufacturer	LILLYCOVER Co., Ltd./KOREA
Date of manufacturing	2022.11
Safety Confirmation certification	ZU10878-18001C
Product composition	Muilli 2(LC-200), C-TYPE Cable /50cm, Manual
KC Safety Certification	R-C-LCV-LC-200

Important Note:

FCC RF Radiation Exposure Statement:

This equipment complies with FCC RF Radiation exposure limits set forth for an uncontrolled environment.


Muilli
by SMARTSKN




Muilli

by SmartSKN

AI Dermoscope

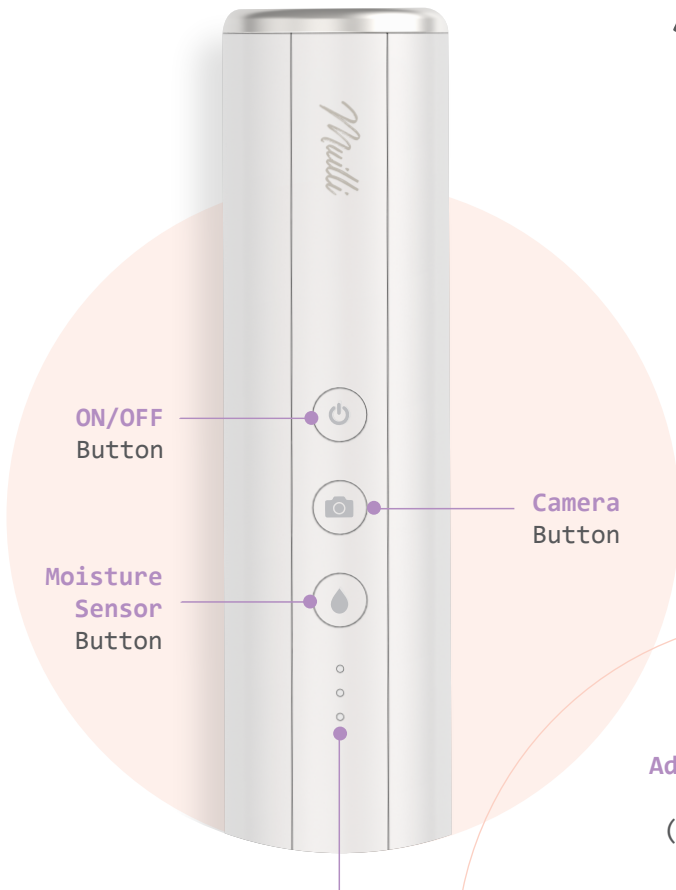

Battery Level Indicator

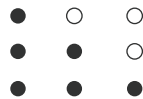

100% 60% 30%

Connect the device to a power source using a USB-C charging cable.

Charging starts when the LED indicators light up.

Advanced Skin Imagery Camera Lens (60x Magnification)



App store
Available on the
App Store



Google Play store
Available on
Google Play



AI Dermoscope Quickstart Guide



Download

Download the smartSKN App on your smartphone or tablet.



App store



Available on the
App Store



Google Play store



ANDROID APP ON
Google Play



Connect

Connect to your smartphone or tablet.



Power:

Press and hold the **ON/OFF button**.

QR Connection:

Scan the **QR code on the bottom** of your Muilli AI Dermoscope device to connect to your device.

Manual Connection:

Use the Settings tab in the app and follow the prompts to connect your device manually.



Setup

Once connected, follow prompts in the app to login or create your account.

See reverse for further instructions



Prep

For the most accurate analysis, remove all the make-up and wait 10-15 mins after washing your face to begin.

Ready?



Scan in 3 Simple Steps



Advanced Skin Imagery
(60x magnification)



Bioimpedance Sensor
(measures moisture-oil level)



Skin Health Questionnaire
(insights into your skin condition and needs)



Results

Real-time Skin Analysis

Receive a diagnostic report of your skin's sensitivity, hydration, aging, imbalances, breakouts, pore size, oil production, pigmentation, and skin barrier health.



Thank you for choosing Muilli AI Dermascope.
For updates, tips, and support, follow us on social media or visit our website at www.smartskn.com.

Need help? Contact our customer support team at support@smartskn.com.

Muilli
by smartSKN